

Product Data

Product Name: Nelaton Catheter


Model: F6, F8, F10, F12, F14, F16, F18, F20, F22, F24

Products Description

Disposable Urethral Catheter is used by patient to drain out the urine. Firstly choose suitable size catheter and lubricate the end with fluid paraffin. Nip the catheter with sterilized hemostats clamp and insert it slowly into the evictions entrance finally to the bladder, then drain out urine.

Intended use: For drainage of urine

Photo of product

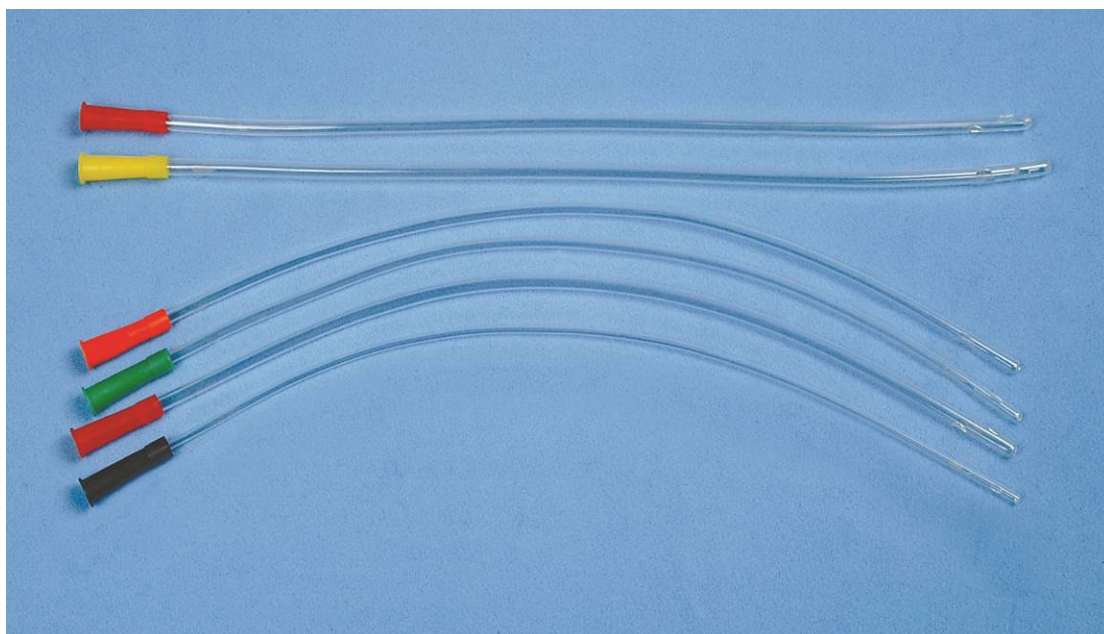

Product configuration

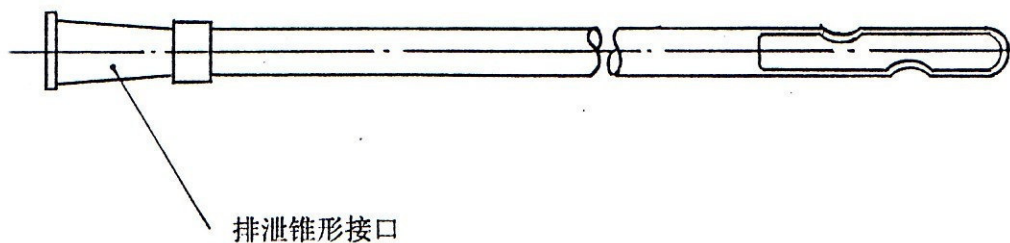


Fig Catheter model illustration (Catheter is made up of tube and connector.)

Product Specification

Specification and basic dimension			Unit mm	
Size	Drainage taper connector color	Whole L	Outer diameter and tolerance	Minimal inner diameter
F6	Aqua	400±10	2.0±0.3	1.05
F8	Cambridge blue		2.7±0.3	1.50
F10	Black		3.3±0.3	2.00
F12	White		4.0±0.3	2.45
F14	Green		4.7±0.3	2.95
F16	Orange		5.3±0.3	3.40
F18	Red		6.0±0.3	3.90
F20	Yellow		6.7±0.3	4.30
F22	Aubergine		7.3±0.3	4.70
F24	Azure		8.0±0.3	5.10